


# Medtronic 5388 dual channel Pacemaker



## Features:

The Medtronic® Model 5388 is a battery-powered, temporary, dual chamber pacemaker designed primarily for temporary antibradycardia pacing therapy. The device provides eight selectable modes of pacing therapy: DDD, DVI, DDI, DOO, VOO, VVI, AOO, and AAI. High-rate, burst pacing therapy up to 800 min<sup>-1</sup> (reciprocal minutes), for atrial tachyarrhythmias, is available in the asynchronous mode.

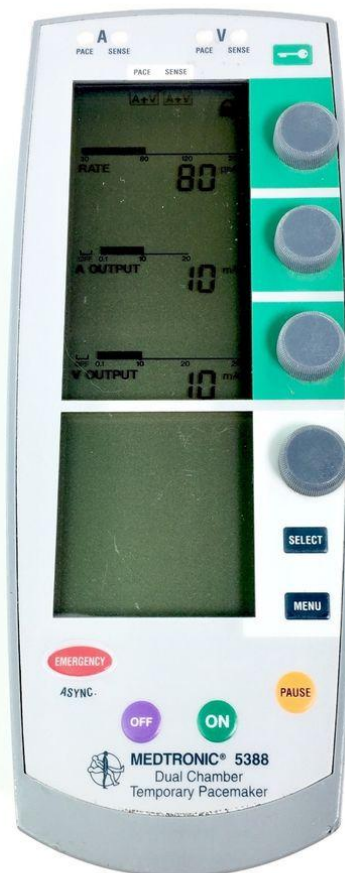

# Features:

- The Medtronic Model 5388 is designed to be reliable, easy to operate, and comfortable to hold.
- Cables with no exposed electrical connections
- Self-test function
- Low Battery indicator
- Lock feature to prevent accidental change of parameters
- Safe, two-step operation to turn the device off
- Runaway rate protection
- Protection from defibrillation shock
- Continuous operation during battery replacement
- Reversible battery polarity
- Electrostatic protection
- Minimized susceptibility to electromagnetic and magnetic interference
- Rubber seals to cover the connector pin receptacles.

## **Medinnova Systems Pvt. Ltd.**

Formerly Soma Tech Pvt. Ltd  
Regd Office: E-114 GIDC – Savli  
Manjusar - Savli Road, At. Po.  
Manjusar - 391 775, Tal. Savli. Dist., Vadodara, Gujarat - India  
Phone: +91-2667-264737  
sales@medinnovasystems.com | www.medinnovasystems.com



---

Note: The technical data given in this publication is for general information and are subject to change without notice. Actual configuration on the unit may vary. Contact our sales representatives for a complete list of details.